



# DosiCare

Operational Dosimetry Software  
Dedicated to Medical Environment



Nuclear  
Power



Healthcare



Homeland  
Security  
& Defense



Labs and  
Education



Industrial and  
Manufacturing

## OVERVIEW

The DosiCare Software easily allows the operational dosimetry management of medical specialists and staff of small and medium-sized infrastructures.

Now, in compliance with new regulations, DosiCare transfers individual dosimetric data, makes available a secured Web interface for consultation of individual dosimetric summaries and makes possible the importation of dosimetric data from an Excel or .csv file.

Thanks to its user-friendly and intuitive interface, DosiCare offers a better follow-up of the dosimetric activity of the agents with a real-time control of their authorizations, the implementation of automated alerts in case of regulatory thresholds being exceeded, the follow-up of the cumulated doses and data transfer.

## KEY FEATURES

- Compliant with new regulations
- User-friendly and intuitive interface
- Web access to essential information (dosimetric data follow-up)
- Automatic EMAIL capabilities
- Complying with regulations
- Integrated database
- Windows compatible
- Consultation of worker data through the Web
- Standard LAPTOP or desktop PC platform
- Possibility of remote network reader management

## ASSOCIATED PRODUCTS


- Readers: LDM 320D, LDM 320W LDM 2000, LDM 3000M
- Dosimeters: DMC 2000S, DMC 2000 GN, DMC 3000, DMC 3000 neutron, DCM 3000 beta

## FUNCTIONAL CHARACTERISTICS

- Consultation of dosimetric summary of an agent via the Web, with secured access
- User-friendly and intuitive interface
- Easy to use (secured Web access to essential information)
- Automatic email to RP person with doses and corresponding data extraction
- Automatic loading of alarm thresholds upon activation
- Graphic display by the worker of his last 30 days upon deactivation
- Customization of the dose and dose rate alarms according to the position occupied or task code used.
- Real-time check of dosimeter activation authorization (medical security clearance, dose limit exceeding)
- Temporary or permanent dosimeter assignment
- Schedule board for:
  - automatic alert for any deactivation in the presence of a dose or dose rate alarm or prealarm
  - automatic alert when a worker comes close to one of his dose limits (configurable)
  - reader connection/link loss
  - active dosimeter followup
- Individual dosimetry management:
  - cumulative totals over the following periods: 1 month, 3 months, 1 year (calendar), 365 days (sliding), 5 years
  - daily dose history
- Report printing:
  - Summaries: individuals, groups of individuals, over a selectable period
  - Collective dosimetry summary per work code
  - Data extraction capability from an Excel file or .csv.

## SYSTEM CONFIGURATION


- Operating system to Windows® 8, on PC and Windows® 2008 Server
- Possible virtualization
- Microsoft SQL Server Express® integrated database
- Basic configuration:
  - 1 PC with DB and dosimeter reader
  - 1 PC for consultation purposes
- Several reader management possibility
- Web consultation station

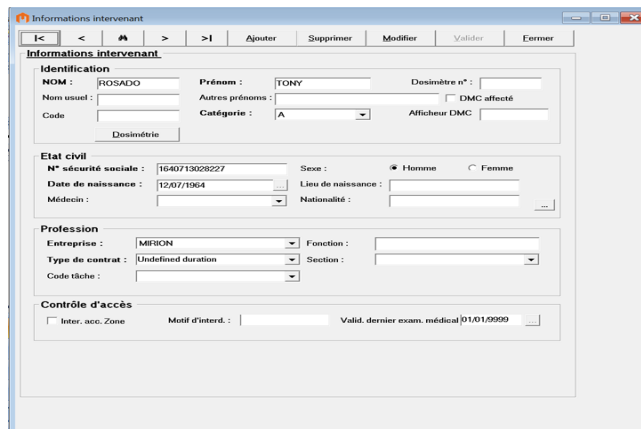


Dosimetric Sheet for a Worker



Screenshot of an entry in controlled area

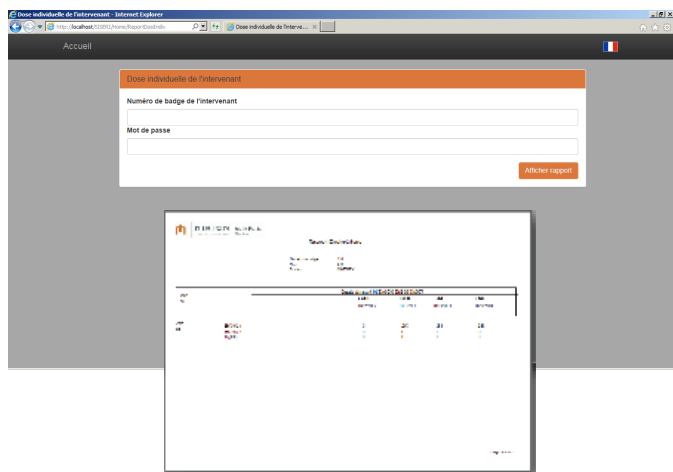

Web consultation screen

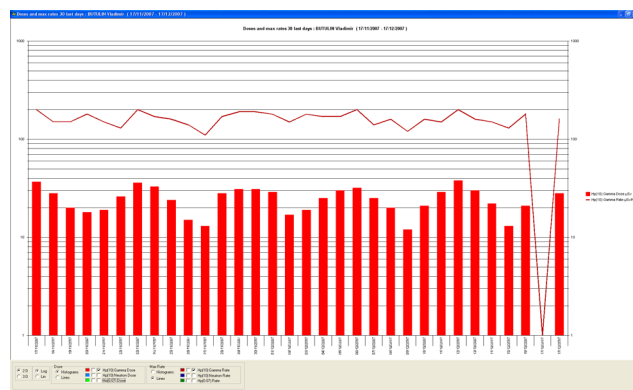

Dosimeter return: statistic display

> CHINA - SHANGHAI  
T: +86 21 6180 6920 | E: info-cn@mirion.com

> FINLAND - TURKU  
T: +358 2 4684 600 | E: info-fi@mirion.com

> FRANCE - LAMANON  
T: +33 (0) 90 59 59 59 | E: info-fr@mirion.com

> GERMANY - HAMBURG  
T: +49 40 85193 0 | E: info-de@mirion.com

> USA - SMYRNA, GEORGIA  
T: +1 770 432 2744 | E: info-us@mirion.com

Copyright (c) 2016 Mirion Technologies, Inc. or its affiliates. All rights reserved. Mirion, the Mirion logo, and other trade names of Mirion products listed herein are registered trademarks or trademarks of Mirion Technologies, Inc. or its affiliates in the United States and other countries. Third party trademarks mentioned are the property of their respective owners.