

ECONO Ultrasound Table

Economically priced, superior ergonomics



Features

- Ergonomic design helps prevent musculoskeletal injuries – comfortable shoulder, arm and hand positions when scanning
- 500 lb patient capacity, suitable for bariatric patients
- Open chassis design for unimpeded patient access
- Drop-down leg section for chair-like comfort
- Height range down to 23" accommodates wheelchair transfers and those with ambulatory difficulties
- Fowler positioning, 0 to 80 degrees, redesigned to accommodate natural body extension and sitting position without slide or shift
- Motorized height and Fowler positioning
- Naugahyde antimicrobial mattress cover with Advanced Beauty Gard top coat provides protection against bacteria
- Durable hand controller activates table positions
- Allows seated or standing work positions

Optional

- Stirrups are designed to cup at the heel for added comfort
- Articulating Scanning Arm Board, adjustable from 0 to 130 degrees
- Headrest for easier thyroid and carotid access
- Folding side rails mount flush to table
- IV Pole easily attaches for quick convenience
- Paper Dispenser accommodates a roll of hygienic table paper. Includes a cutter strap.

Description


The Econo Ultrasound Table is designed on the same platform as the Ultra Pro™ Table with similar features and functions. A more cost-effective alternative, it is offered without Trendelenburg positioning or central floor-locking system.

The same ergonomics have been incorporated into the design of the Econo Ultrasound Table, making it easy for a sonographer to achieve quality images without discomfort to the patient or sonographer.

Motorized height and Fowler positioning which accommodates natural body extension and sitting position, a height range as low as 23" to allow wheelchair transfers, and a hand controller are just some of the features that make the imaging experience easy with the Econo Ultrasound Table.

The Econo Table can be customized to suit your department needs with an assortment of table options. Flush-mounted side rails enable closer patient access without sonographer contortion. An articulating scanning arm board adjustable from 0 to 130 degrees, and a headrest that drops down 15 degrees to allow for carotid studies are just some of the items which can be added to the table.

Biodex | Econo Ultrasound Table

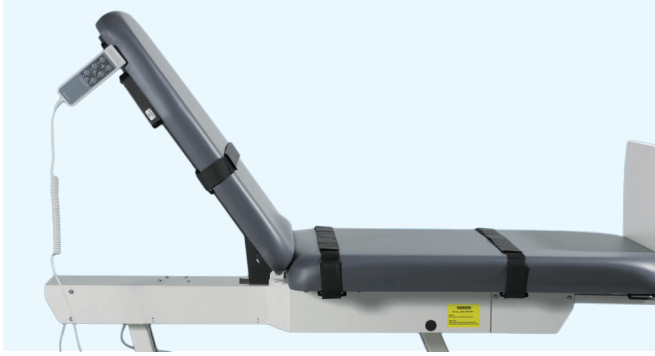

The motorized Fowler back has minimal patient shift of no more than 2-3 inches.

Specifications

• Dimensions:

Overall: 70" l x 30" w (177.8 x 76.2 cm); 35" w (88.9 cm) with optional side rails

• **Weight:** 335 lb (152 kg)

• **Drop-Down Leg Section:** 12.6" l (32 cm)

• Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Fowler Positioning: 0° to 80° infinitely adjustable

• Controls:

Hand Control: Activates height and Fowler back

• **Tabletop:** Three primary sections: torso section with Fowler positioning, center section remains fixed and leg section drops down to 40° and 80° for stirrup access

• **Wheels:** 5" (12.7 cm) Individual locking swivel casters

• **Patient Restraints:** Two body straps

• **Mattress:** Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; Naugahyde antimicrobial mattress cover with Advanced Beauty Gard top coat provides protection against bacteria

• **Upholstery Color:** Graphite

• **Finish:** Powder coat

• **Patient Capacity:** 500 lb (227 kg); weight tested to four times the patient load rating

• **Power:** 115 VAC or 230 VAC

• **Warranty:** Three years parts and labor

• **Certifications:** SGS listed to UL60601-1 and CAN/CSA C22.2 No.601.1.M90 and EN 60601-1 standards

• All Ultrasound Tables are compatible with Hoyer Lifts

• **Primary Function:** General Ultrasound



Optional flush-mounted side rails provide added safety without obstructing patient access.

058-726 Table, Ultrasound, Econo, 115 VAC
Includes hand controller.

058-727 Table, Ultrasound, Econo, 230 VAC
Includes hand controller.

Related:

058-704 Chair, Ergonomic, Sonography

Optional:

058-633 Side Rails, Folding

058-652 Stirrups, Retractable

058-736 Arm Board, Articulating Scanning

058-737 IV Pole

058-611 Dispenser, Paper

058-738 Headrest

Ordering Information



mirion.com/biodex