

ECHO/VASC PRO™

Vascular Echocardiography Table

Achieve quality vascular and echocardiography images



Features

- Ergonomic design helps prevent musculoskeletal injuries – comfortable shoulder, arm and hand positions when scanning
- Extra-wide top with 500 lb patient capacity, suitable for bariatric patients
- Open chassis design for unimpeded patient access
- Cardiac drop-down cushion, allows release and return from either side of table
- Industry's largest cardiac scanning cutout – 15.75" l x 9.5" w
- Reverse Trendelenburg to 30 degrees
- Height range down to 23" accommodates wheelchair transfers and those with ambulatory difficulties
- Fowler positioning, 0 to 80 degrees, redesigned to accommodate natural body extension and sitting position without slide or shift
- Motorized height, Fowler positioning, Trendelenburg motions and auto-level adjustment
- Adjustable foot support with comfortable access to calves
- Naugahyde antimicrobial mattress cover with Advanced Beauty Gard top coat provides protection against bacteria
- Central floor-locking system maintains table stability
- Durable hand controller activates table positions
- Allows seated or standing work positions

Optional

- Foot Controller is ideal for when hands are busy with other functions
- Folding side rails mount flush to table
- Articulating Scanning Arm Board, adjustable from 0 to 130 degrees
- IV Pole easily attaches for quick convenience
- Paper Dispenser accommodates a roll of hygienic table paper. Includes a cutter strap.
- Headrest drops down 15 degrees to allow for carotid and thyroid studies

Description

The Echo/Vasc Pro Vascular Echocardiography Table is the logical choice for performing vascular scanning and echocardiography. It is designed with many features that help obtain quality images and keep the sonographer safe from injury.

The Echo/Vasc Pro Table features 30 degree reverse Trendelenburg positioning, ideal for expanding lower extremity veins for maximum vessel dilation. The adjustable foot support provides the sonographer with comfortable access to patients of various heights, without crouching. When not in use, the foot support stores beneath the table, allowing clear access to the tabletop. Body straps stabilize the patient and ensure safety throughout the entire procedure. The design also includes a cardiac drop-down cushion, with the ability to release and return from either side of the table, providing ample scanning access for an unobstructed apical approach.

As with all Biodex™ Ultrasound Tables, the Echo/Vasc Pro Table features Fowler positioning to accommodate natural body extension and sitting position without slide and shift, optional side rails for easy access to the patient, an optional Articulating Scanning Arm Board adjustable from 0 to 130 degrees, and soft antimicrobial vinyl upholstery.

Biodex | Echo/Vasc Pro Vascular Echocardiography Table



Adjustable every 2", up to 14", the foot support provides the sonographer comfortable access to patients of various heights.

Specifications

• Dimensions:

Overall: 73" l x 30" w (185.4 x 76.2 cm); 35" w (88.9 cm) with optional side rails

• **Weight:** 335 lb (152 kg)

• **Cardiac Scanning Cutout:** 15.75" l x 9.5" w (40 x 24.1 cm) drop-down cushion releases from either side of the table

• **Sonographer Cutout:** 14.75" l x 8" w (37.4 x 20.3 cm) adjustable sonographer cushion drops down or folds up for maximum access to the patient

• Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to +15°, 0° to -30°

Fowler Positioning: 0° to 80° infinitely adjustable

• Controls:

Hand Control: Activates height, Trendelenburg motions, Fowler positioning and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler positioning

• **Tabletop:** Two primary sections: torso section features sonographer and patient cutouts with Fowler positioning; lower body section remains fixed

• **Wheels:** 5" (12.7 cm) Individual locking swivel casters; central floor-locking system


• **Patient Restraints:** Two body straps

• **Foot Support:** Adjustable every 2", up to 14". Stores beneath the table when not in use

• **Mattress:** Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; Naugahyde antimicrobial mattress cover with Advanced Beauty Gard top coat provides protection against bacteria

• **Upholstery Color:** Graphite

• **Finish:** Powder coat

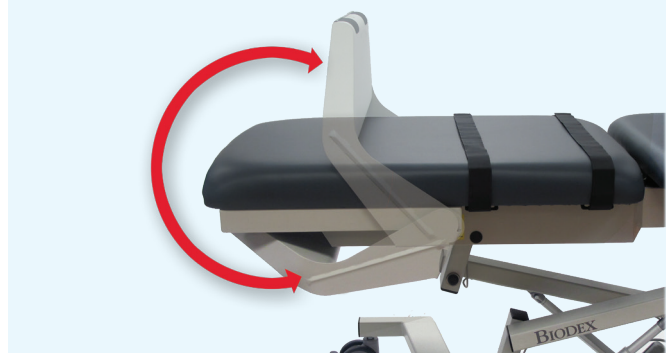

When not in use, the foot support stores beneath the table allowing clear access to the tabletop.

- **Patient Capacity:** 500 lb (227 kg); weight tested to four times the patient load rating
- **Power:** 115 VAC or 230 VAC
- **Warranty:** Three years parts and labor
- **Certifications:** SGS listed to UL60601-1 and CAN/CSA C22.2 No.601.1.M90 and EN 60601-1 standards
- All Ultrasound Tables are compatible with Hoyer Lifts
- **Primary Function:** Echocardiography, Vascular

Ordering Information

- | | |
|---------|--|
| 058-702 | Table, Vascular Echocardiography, Echo/Vasc Pro, 115 VAC
Includes hand controller and foot support. |
| 058-707 | Table, Vascular Echocardiography, Echo/Vasc Pro, 230 VAC
Includes hand controller and foot support. |

Related:

- | | |
|---------|------------------------------|
| 058-704 | Chair, Ergonomic, Sonography |
|---------|------------------------------|

Optional:

- | | |
|---------|----------------------------------|
| 058-741 | Foot Controller |
| 058-633 | Side Rails, Folding |
| 058-736 | Arm Board, Articulating Scanning |
| 058-737 | IV Pole |
| 058-611 | Dispenser, Paper |
| 058-738 | Headrest |



mirion.com/biodex