




ConnectStudio™


Connectivity Management Software

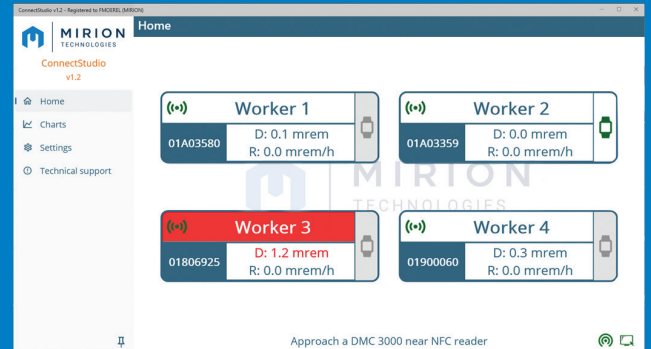

ConnectStudio software strengthens exposure monitoring and supervision by enabling connectivity of the RDS-Med™ and DMC 3000™ dosimeter with Bluetooth® module, and related ecosystem.

RELATED PRODUCTS

- DMC 3000 with Bluetooth Module
- RDS-Med Survey Meter and Alarm Box

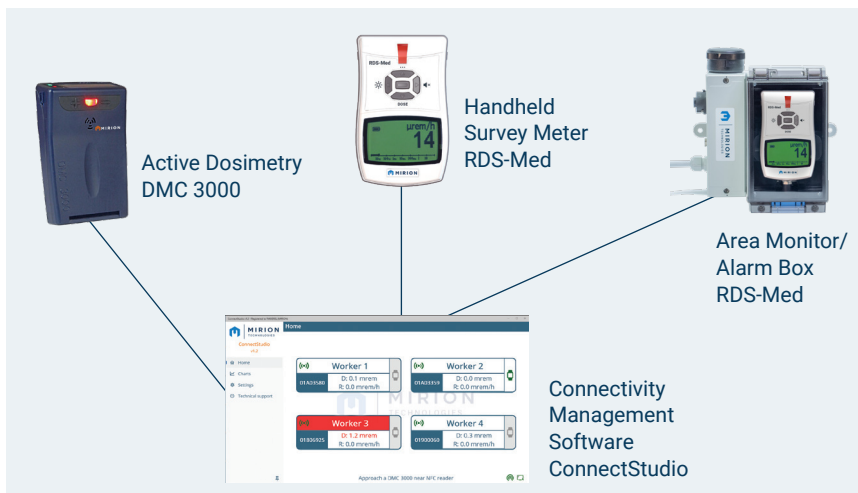

FEATURES


- ✓ Efficient and easy pairing management
- ✓ Operational telemetry supervision at a glance
- ✓ Visual display of dose and dose rate to support dose exposure reduction for procedures or processes
- ✓ Credit dose display to support dose exposure reduction for procedures or processes
- ✓ Facilitates wireless, short-range supervision
- ✓ Compatible WiFi network
- ✓ Segment work into areas such as radiopharmacy, hot labs, radiology and emergency department



FUNCTIONALITY

To simplify the pairing and connection with Mirion connected objects

ConnectStudio software allows for an easy pairing with the NFC reader, and to the DMC 3000 with Bluetooth Module and RDS-Med Survey Meter. The software manages the connections, and reconnections automatically in case of a lost signal.

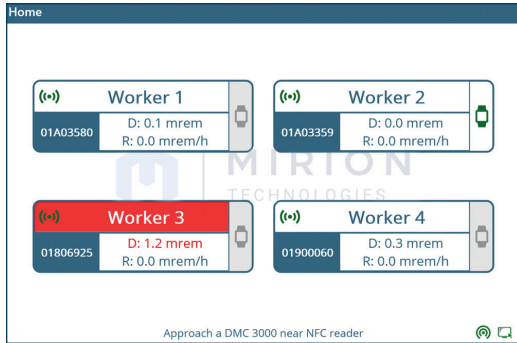

Gateway between Bluetooth connection and Mirion WRM network

ConnectStudio software provides a bridge to the legacy Mirion WRM system (remote supervision) through WiFi network.

RECOMMENDED & MINIMUM COMPUTER CONFIGURATION


- Tablet or PC - Windows 10, version 1809 minimum
- Android version 12
- 10" display
- 8 GB RAM
- Bluetooth Low Energy 4.2

The maximum number of simultaneously connected dosimeters to ConnectStudio software can vary from six to more than 15 depending on tablet limitations.

- Standard tablet:
 - HP Pro x2, HP Elite X2, Microsoft Surface Pro 7
- NFC reader:
 - Any NFC reader integrated on a compatible tablet
 - Any NFC reader configurable as a smartcard reader - ACS ACR122U; ACS ACR1252U

OPERATIONAL TELEMETRY SUPERVISION

ConnectStudio software offers a quick operational display of data from connected products. The information can be displayed with a graphical view, accumulated dose view (dosimeter) or tile view. This functioning view allows for any necessary quick action in case of an alarm or pre-alarm (dose/dose rate/time alarm).

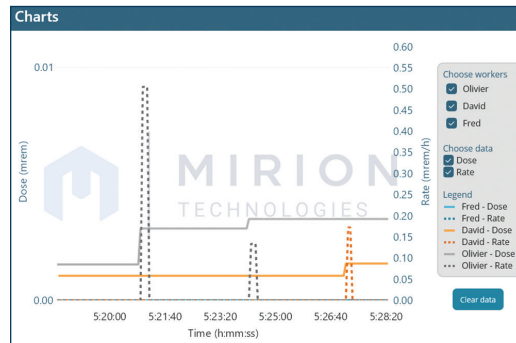

COMPATIBILITY

- DMC 3000 Bluetooth Low Energy
- RDS-Med Survey Meter and probes
- RDS-Med AlarmBox

Available Through Mirion Medical Companies: Capintec & Sun Nuclear

Health Physics for Healthcare solutions enable radiation monitoring in healthcare applications — a critical aspect of safety and compliance for all facilities using radiation in patient diagnosis and care.

